Posey® Twice-As-Tough® Cuffs for Stretcher 2750, 2755 Application Instructions for Wrist and Ankle

Rx ONLY

DESCRIPTION OF PRODUCT: Padded limb restraints with double-security closure. For stretcher application only.

INTENDED USE:

- Patients assessed to be at risk of disrupting life-saving treatments (e.g., chronic tube pulling) or in danger of injury to themselves or to others.
- Follow your hospital's restraint policies and procedures which are in compliance with USA CMS guidelines and state laws, or other governing agencies outside the USA.


CONTRAINDICATIONS:

Do not use this device with someone who has continued highly aggressive or combative behavior, self-destructive behavior, or deemed to be an immediate risk to others or to self. Clinicians may need to use additional interventions in conjunction with restraints.

APPLICATION INSTRUCTIONS:

Follow these steps to apply device (repeat steps 1-4 on each side):


1. Place restraint on the stretcher at both the wrist position and ankle position of the stretcher.
2. Bring the ends of the straps down through the inside of the side rails so they do not interfere with the movement of the side rails. Thread the end of one strap over the top, around the frame, and first attach the black hook-and-loop together then the blue or red hook-and-loop together (Fig. 1). Be sure to overlap at least one inch (3 cm).
3. Tug the restraint to confirm it is secure to the stretcher.
4. Position the patient on the stretcher.

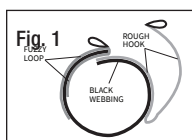APPLYING THE CUFFS:

(Repeat steps 1-2 for each side):

1. Wrap the neoprene piece around the wrist/ ankle so the buckle and connecting strap is on the ulnar side of the wrist or lateral malleolus of the ankle. Attach the black hook-and-loop pieces together and pull back the black pull tab on the blue or red fuzzy piece, followed by the blue or red hook-and-loop pieces. The fuzzy piece should be sandwiched between the two pieces of hook (Fig. 2). Be sure to overlap at least one inch (3 cm).
2. Press the hook-and-loop closure together firmly and make sure it adheres securely. Slide ONE finger (flat) between the cuff and the inside of the patient's wrist/ankle to ensure proper fit.

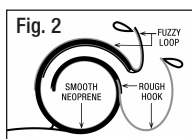POSEY TWICE-AS-TOUGH CUFFS:

REF 2750	Stretcher, Regular/Blue Size
REF 2755	Stretcher, Large/Red Size

PRECAUTIONS:

- A clinical assessment and decision are required when used with monitoring lines or if patient has a wound or dislocated/fractured limb.
- Check the patient regularly to ensure that circulation is not impaired. Serious injury may occur if the cuffs restricts circulation when the limb holder is applied.
- Before each use, check cuffs and straps for cracks, tears, and/or excessive wear or stretch, broken buckles or locks, and/or that hook-and-loop adheres securely as these may allow patient to remove cuff. Discard if device is damaged or if unable to lock.
- **WARNING:** Additional or different body or limb restraints may be needed (visit tidiproducts.com):
 - If the patient pulls violently against the bed straps.
 - To reduce the risk of the patient getting access to the line/wound/tube site.
 - To prevent the patient from flailing or bucking up and down and causing self-injury.

BED SAFETY:

Refer to the Food and Drug Administration (FDA), or other governing agencies outside the USA, for the most recent Hospital Bed Safety Guidelines and the bed manufacturer's instructions for use.

EN

Report any serious incident related to device use to TIDI Products and the Member State Competent Authority. | Translations of this document in French, Spanish, Dutch, German, Italian, and Portuguese are located at www.tidiproducts.com/IFU.

ADDITIONAL SAFETY AND LAUNDERING INSTRUCTIONS ON OTHER SIDE.

TIDI Products, LLC • 570 Enterprise Drive, Neenah, WI 54956 USA
Phone: 1.800.447.6739 • International: +1.920.751.4036 • www.tidiproducts.com

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UK Responsible Person
Emergo Consulting (UK) Limited
c/o CR360 – UL International
Compass House, Vision Park Histon
Cambridge CB249BZ
United Kingdom

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UK CA UK Persona Responsible: Emergo Consulting (UK) Limited c/o CR360 – UL International Compass House, Vision Park Histon Cambridge CB249BZ Reino Unido

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
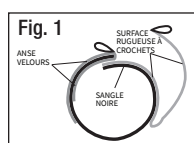
APPLICATION DES ATTACHE-POIGNETS/CHEVILLES:

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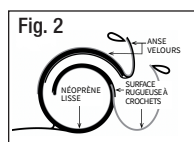ATTACHE-POIGNETS/CHEVILLES POSEY TWICE-AS-TOUGH:

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MDSS GmbH
Schiffgraben 41
30175 Hanovre,
Allemagne

UK Responsible
Emergo Consulting (UK) Limited
c/o CR360 – UL International
Compass House, Vision Park Histon
Cambridge CB249BZ
Royaume-Uni

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UK Pessoa responsável
Emergo Consulting (UK) Limited
c/o CR360 – UL International
Compass House, Vision Park Histon
Cambridge CB249BZ
Reino Unido

19232 REV F 2021-08-02


For RX PRODUCTS ONLY. Federal law (USA) restricts RX devices to sale by or on order of a physician. For use in a licensed healthcare facility only.

STAFF TRAINING: TIDI offers in-service training aids at no charge. Contact TIDI online at www.tidiproducs.com or call toll-free at 1.800.447.6739.

BEFORE APPLYING ANY RESTRAINT: Follow your hospital's restraint policies and procedures that are compliant with USA CMS guidelines and state laws or other governing agencies outside the USA.

⚠️ ADDITIONAL WARNINGS:

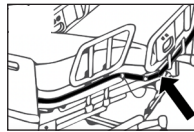1. ALWAYS monitor patient per facility policy. Improper application or use of any restraint may result in serious injury or death. Be aware that constant monitoring may be required for:
 - Aggressive or agitated patients; and
 - Patients in a supine position; ensure that being in this position would not cause any complications for those who are at risk of vomiting while being restrained.
2. NEVER alter or repair this product. ALWAYS inspect before each use: Check for broken stitches or parts; torn, cut or frayed material; or locks, buckles, or hook-and-loop fasteners that do not hold securely. DO NOT use soiled or damaged products. Doing so may result in serious injury or death. Dispose of damaged products per facility policy for BIOHAZARDOUS material.
3. ALWAYS secure straps to a part of the frame that moves with the patient (not foot/head board or siderails), out of the patients reach, using quick-release ties (see drawing below) or buckles. These allow easy release in the event of an accident or fire.
 

Test to make sure straps cannot tighten, loosen, or slip and create excess slack. If this occurs, the patient may slide off the chair or bed, increasing the risk of serious injury or suffocation. Restraint release is an important part of facility fire and disaster drills. Straps can be cut with scissors in an emergency.
4. NEVER secure restraint strap to side rail or head/foot board.
5. NEVER use Posey products on toilets, or on any chair or furniture that does not allow proper application as directed in the Application Instructions. DO NOT use at home.
6. NEVER expose this product to open flame, fire, smoking materials, or high heat sources. Some products may melt or ignite and burn. The facility smoking/no smoking policy should be strictly enforced.
7. NEVER use a Posey product as a seat belt in a moving vehicle. Posey products are not designed to withstand the force of a crash or sudden stop.
8. Test Zippers or hook-and-loop fasteners before each use. DISCARD device if it does not fasten securely.

LAUNDERING INSTRUCTIONS (if applicable):

- Fasten all buckles and locks to reduce risk of damage during wash and dry cycles. DO NOT put buckles or locks through extractors. For maximum life, launder in a laundry bag.
- Stainless Steel Locks: Posey recommends that locks be treated with a dryfilm lubricant (silicone) after each cleaning. For best results, spray a small amount into the lock. Work the action of the lock with a key several times to ensure that all the surfaces are well lubricated.
- Before laundering, zip up and turn the product inside out to protect zipper.
- Hook-and-loop fasteners may collect lint after repeated use or laundering, reducing grip strength. Fasten the "hook" to the "loop" before laundering to help prevent lint buildup. As needed, use a stiff-bristle brush to remove lint from the "hook" side.

• These products, other than foam products, can be machine washed under CDC guidelines for material soiled with blood or bodily fluid.

- For non-contaminated material, use lower temperature wash and dry cycles to extend product life.

• For foam products:

- Make sure products are completely dry before use.

STORAGE AND HANDLING:

- This device is designed for use in normal indoor environments.
- This device may be stored in ambient warehouse temperatures at normal humidity levels. Avoid excess moisture or high humidity that may damage product materials.

HOW TO TIE THE POSEY QUICK-RELEASE TIE

1. Wrap strap once around the frame that moves with the patient (do not attach to side rail or head/footboard), leaving at least an 8" (20 cm) tail. Fold the loose end in half to create a loop and cross it over the other end.
2. Insert the folded strap where the straps cross over each other, as if tying a shoelace. Pull on the loop to tighten.
3. Fold the loose end in half to create a second loop.
4. Insert the second loop into the first loop.
5. Pull on the loop to tighten. Test to make sure strap is secure and will not slide in any direction.
6. Repeat on other side. Practice quick-release ties to ensure the knot releases with one pull on the loose end of the strap.

EN

Report any serious incident related to device use to TIDI Products and the Member State Competent Authority. | Translations of this document in French, Spanish, Dutch, German, Italian, and Portuguese are located at www.tidiproducs.com/IFU.

DISPOSAL: Dispose of product according to local or facility guidelines if product is worn or damaged.

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© Posey Products, LLC. Tous droits réservés.
TIDI Products, LLC • 570 Enterprise Drive, Neenah, WI 54956 EUA